



Advanced Combinational Clinical Therapy Rack

Designed to support modern treatment workflows, recovery-focused care, mobility support, and expanded non-invasive therapy offerings.



RSWT

Radial
Shockwave
Therapy



FSWT

Focus
Shockwave



VIBE

Vibration
Therapy



ULTRA

Ultrasound
Therapy



PMST

Magneto
PMST



HILT

High Intensity
Laser Therapy



LLLT

Low-Level
Laser Therapy



IFC

Electrotherapy

Content

What is NerveCube?	1
PMST	5
Electromagnetic Recovery & Neuromuscular Support	5
VIBE	5
Circulation, Recovery & Muscle Stimulation	5
FOCUS SHOCKWAVE	6
Deep Tissue Recovery & Mobility Support	6
RSWT	6
Soft Tissue Recovery & Rehabilitation Support	6
Ultrasound Therapy	7
Targeted Tissue Stimulation & Recovery Support	7
HILT	7
High-Intensity Therapy for Recovery Applications	7
LLLT	8
Cellular Stimulation & Soft Tissue Support	8
TECAR Therapy	8
PMST LOOP	9
Common Clinical Applications	10
User Interface	10
Therapy Selection	10
Treatment Control	10
Protocol Management	10
Detailed Specifications	14
Available by Special Order	14
Detailed Specifications	14

What is NerveCube?

NerveCube is a combinational therapy platform that integrates advanced treatment technologies into a single clinical system. Designed for chiropractors, rehabilitation clinics, sports recovery providers, and modern therapy centers, NerveCube enables flexible treatment combinations through advanced therapy modules.

The platform combines Magneto (PMST), Shockwave, High-Intensity Laser, Ultrasound, Vibration, and Cold Laser technologies to support adaptable clinical workflows across a wide range of rehabilitation, recovery, mobility-focused, neuromuscular, and performance-support applications. By consolidating complementary modalities into a unified architecture, NerveCube helps clinics streamline workflow efficiency, expand treatment versatility, improve the patient experience, and reduce the need for multiple standalone systems.

Advanced Combo 1

			
RSWT	FSWT	VIBE	ULTRA
Radial Shockwave Therapy	Focus Shockwave	Vibration Therapy	Ultrasound Therapy

Advanced Combo 2

			
PMST	HILT	LLLT	IFC
Magento PMST	High Intensity Laser Therapy	Low-Level Laser Therapy	Electrotherapy





Clinical Grade Combinational Treatment Stack



Ten Advanced Technologies




VIBE



FOCUS



PULSE

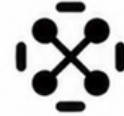

IFC



PMST



ULTRA



HILT



LLLT

ALSO AVAILABLE



TECAR



LOOP



PMST

Magneto PMST

Electromagnetic Recovery & Neuromuscular Support

PMST systems generate controlled pulsed electromagnetic fields designed to support tissue stimulation, circulation, recovery-focused care, and rehabilitation applications.

These electromagnetic fields may help support neuromuscular stimulation, mobility-focused care, nerve-related rehabilitation applications, and recovery-oriented treatment protocols across a wide range of clinical environments.




Tesla 7



Tesla 5

Nerve Cube

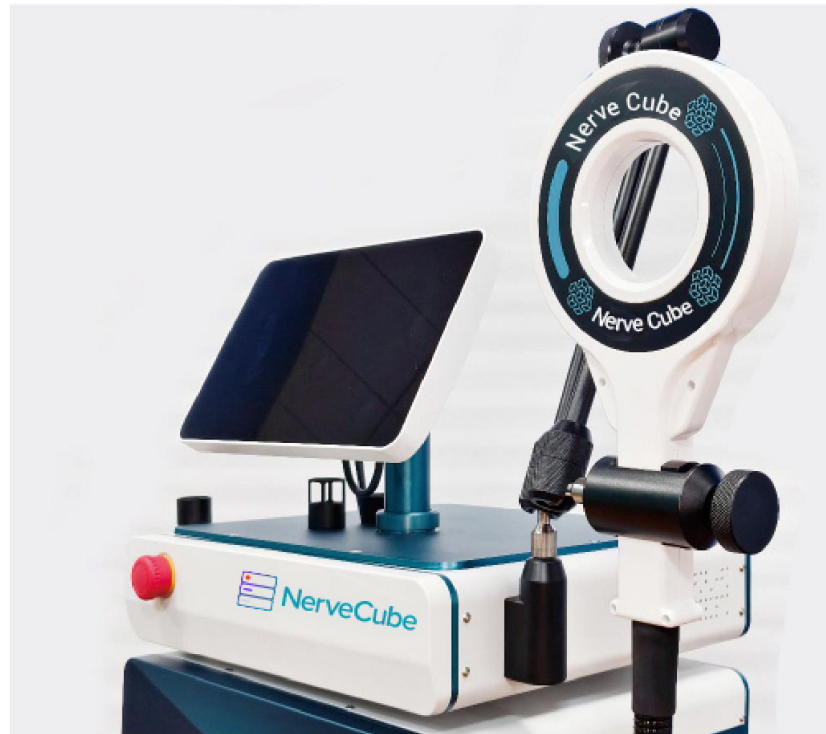

VIBE

Vibration Therapy

Circulation, Recovery & Muscle Stimulation

Vibration therapy delivers high-frequency mechanical oscillations through handheld or platform-based devices. These controlled vibrations may help support circulation, reduce muscle tension, and assist rehabilitation and recovery protocols.

The mechanical stimulation is commonly utilized in chiropractic, rehabilitation, sports recovery, mobility-focused care, and soft-tissue recovery applications.



FOCUS SHOCKWAVE

FSWT

Deep Tissue Recovery & Mobility Support

Focused shockwave technology utilizes high-voltage electromagnetic coils to generate precise, high-energy acoustic waves.

Unlike radial shockwave systems, focused shockwave delivers concentrated energy to targeted tissue depths with high precision, making it suitable for deep-tissue rehabilitation, recovery-focused care, mobility support, and soft-tissue applications.



Nerve Cube



RSWT

Radial Shockwave Therapy

Soft Tissue Recovery & Rehabilitation Support

RSWT utilizes high-energy pressure waves generated pneumatically through the handpiece applicator.

These acoustic waves are designed to support natural rehabilitation processes and are commonly utilized across chiropractic, orthopedic, sports recovery, mobility-focused care, and rehabilitation environments.

Ultrasound Therapy


Targeted Tissue Stimulation & Recovery Support

Ultrasound therapy uses high-frequency sound waves delivered through a handheld applicator to provide targeted tissue stimulation.

The acoustic energy may help support circulation, rehabilitation protocols, tissue relaxation, joint-focused applications, and recovery-oriented treatment programs.



Nerve Cube

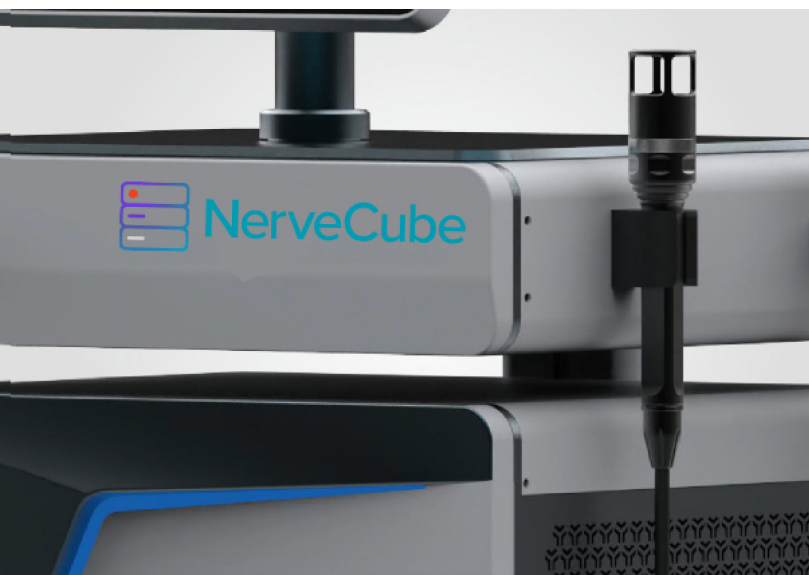

HILT

High-Intensity Laser Therapy

High-Intensity Therapy for Recovery Applications

NerveCube HILT delivers high-intensity laser energy into deep tissues to support rehabilitation protocols involving muscles, tendons, joints, and nerves.

The system is designed to support deep-tissue stimulation, joint-focused rehabilitation, mobility support, recovery-oriented care, and rehabilitation applications.


LLLT

Low-Level Laser Therapy

Cellular Stimulation & Soft Tissue Support

LLLT utilizes 808 nm / 1064 nm laser wavelengths to support cellular stimulation and rehabilitation applications without significant thermal impact.

The modality is commonly used in recovery-focused care, rehabilitation applications, soft-tissue support, targeted treatment protocols, and mobility-oriented programs.

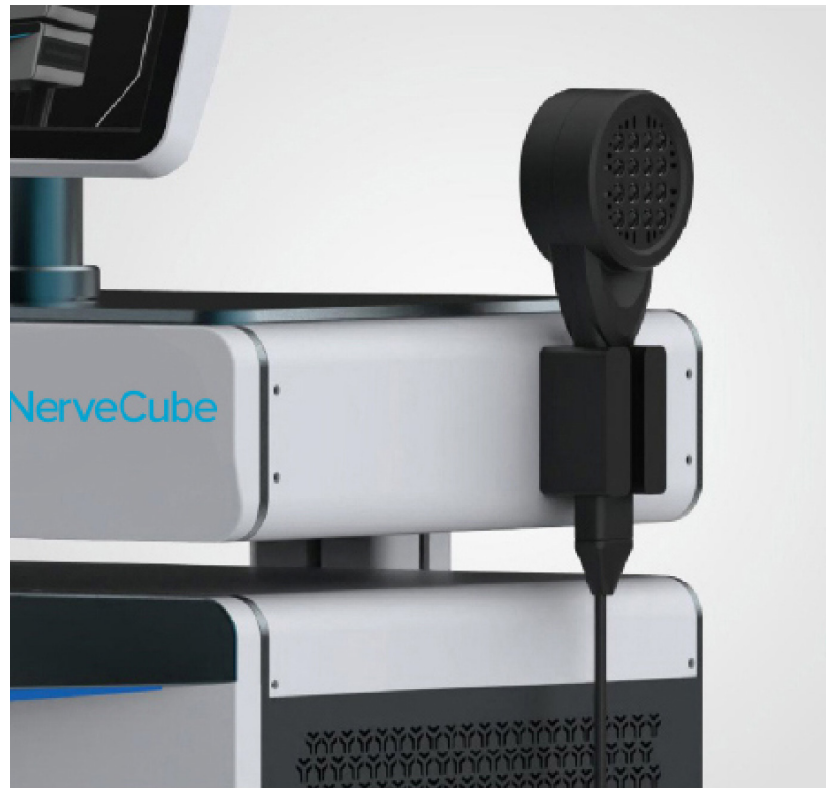

Nerve Cube



TECAR Therapy

Available by special order — Pulse TECAR · PEMF LOOP

TECAR (Capacitive & Resistive) delivers radiofrequency energy through electrodes to generate deep, targeted heat. The RF current stimulates tissue metabolism, improves circulation, and supports recovery in muscles, joints, and ligaments.

PMST LOOP

Available by special order — LOOP

PMST LOOP delivers high-intensity, precisely controlled PEMF pulses for deeper tissue penetration. This non-invasive therapy supports cellular repair, reduces inflammation, and promotes pain relief. It is suitable for both rehabilitation and wellness applications.



Nerve Cube

Common Clinical Applications

- Neuromuscular Stimulation
- Recovery-Focused Care
- Mobility Support
- Joint & Soft Tissue Applications
- Rehabilitation Support
- Circulation Support
- Non-Invasive Pain Management Applications
- Performance Recovery Protocols
- Recovery-Oriented Treatment Workflows
- Mobility-Focused Rehabilitation Programs

User Interface

Therapy Selection

Intuitive touchscreen dashboard providing rapid access to therapy modules, protocols, and treatment configurations.

Treatment Control

Adjust intensity, frequency, pulse settings, and treatment duration for precise therapy customization.

Protocol Management

Save, edit, and manage customized treatment protocols for efficient clinical workflows and operational consistency.



Vibration interface



Radial shockwave interface



Focus shockwave interface



COMBINATIONAL DESIGN

1) NerveCube – Advanced Combo #1


Control Box

- Vibration
- Ultrasound

Therapy Box

- Focused Shockwave
- Radial Shockwave

Supports mobility-focused care, soft tissue recovery, circulation support, neuromuscular stimulation, and rehabilitation-oriented workflows.

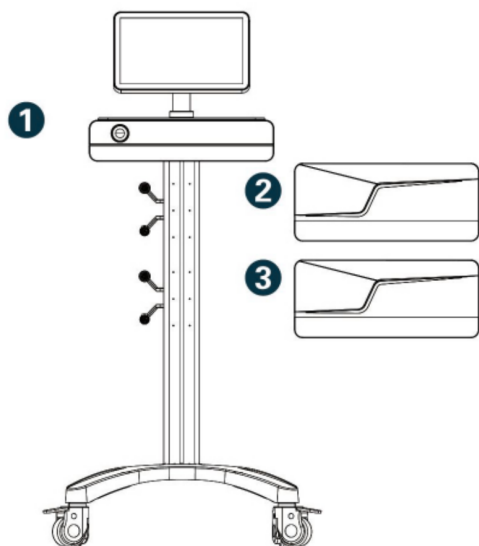

2) NerveCube – Advanced Combo #2


Control Box

- Cold Laser
- Electrotherapy

Therapy Box


- HILT
- Magneto-EMTT

Designed to support recovery-focused applications, tissue stimulation, mobility support, joint-focused rehabilitation, and modern performance recovery workflows.

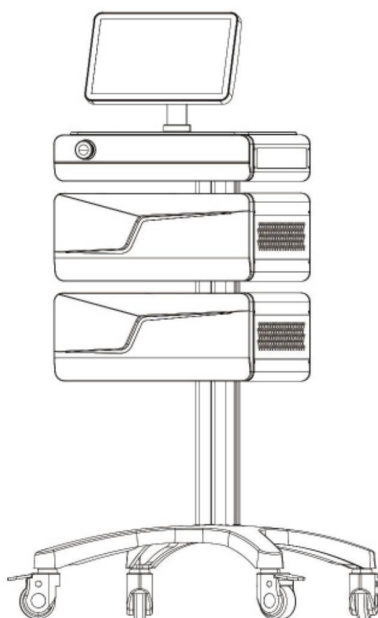


AVAILABLE BY SPECIAL ORDER


- Pulse TECAR
- PEMF LOOP

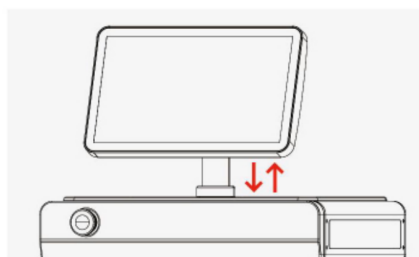

Height-adjustable screen

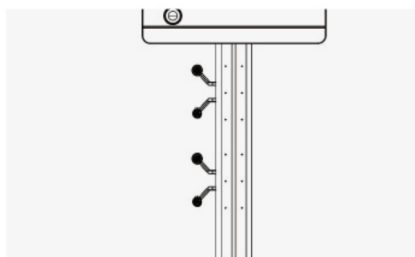

Sturdy support stand

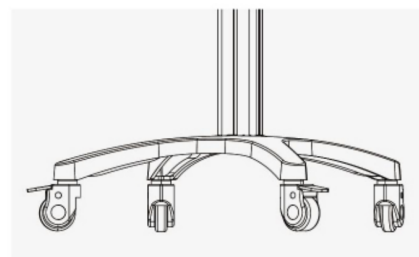

Smooth 360° caster wheels for easy mobility

ASSEMBLY STEPS

1. Install the Stand Assembly

Secure the base and upright support structure.

2. Mount the Therapy Modules

Attach the Control Box and Therapy Box modules.


3. Connect Cables & Accessories

Connect handpieces, applicators, and system accessories.


4. Power On & Configure

Initialize the system and configure treatment settings.


Head Unit & Display Installation (connector attachment)

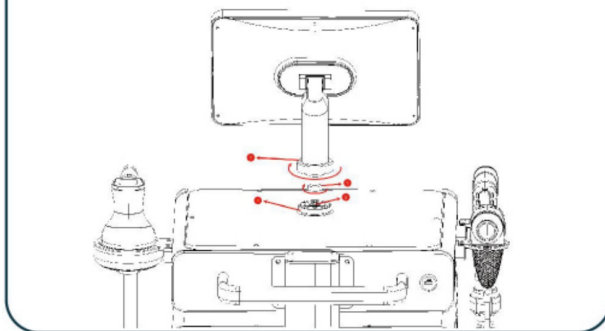




Box Data Connection

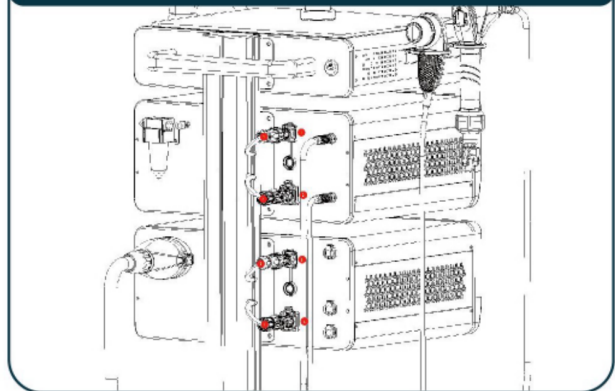

Base Installation (4 screws)

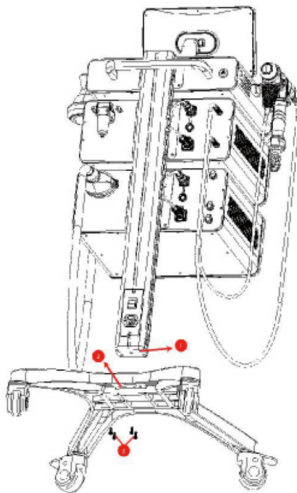

Installation finished

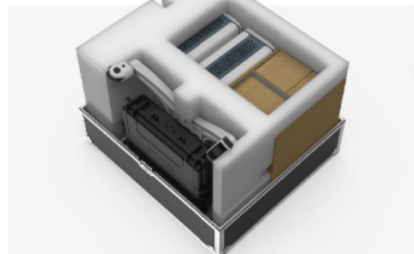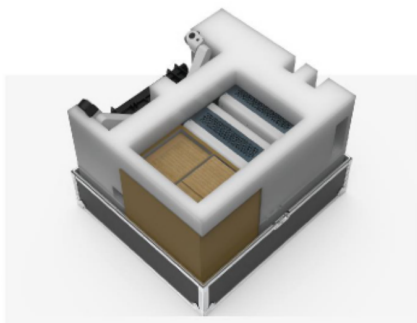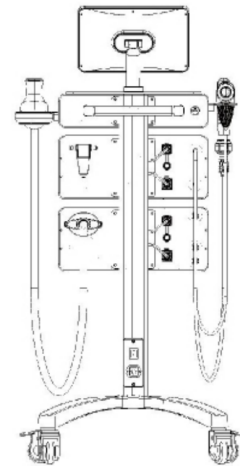

Secure Protection · Damage-Free Delivery

Detailed Specifications

FOCUS SHOCKWAVE	
Intensity	0.03-0.39 mJ/mm²
Frequency	1-10 Hz
Therapeutic depth	120 mm
Pressure	5-50 MPa

MAGNETO	
Intensity	4-7 Tesla
Frequency	1-150 Hz & 1-3 kHz
Therapeutic depth	18 cm
Mode	ST, MT, LT, EXPERT

RADIAL SHOCKWAVE	
Intensity	0.5-10 Bar
Frequency	1-21 Hz

VIBRATION	
Intensity	0.5-2.5 Bar
Frequency	50 Hz

LASER	
Output power	45 W (30 W + 15 W)
Wavelength	980 nm + 1064 nm
Frequency	0.2-20000 Hz
Pulse width	10 µs - 3 s

COLD LASER	
Output power	5600 mW (5.6 W)
Wavelength	650 nm & 808 nm
Mode	Continuous & Pulse

IFC & TENS	
Carrier frequency	1-5 kHz
Modulation frequency	1-100 Hz
Mode	Continuous & Pulse, 0.1-2 Hz

ULTRASOUND THERAPY	
Frequency	1 MHz & 3 MHz
Intensity	0.2-3 W/cm²
Mode	Continuous & Pulse

Available by Special Order

PEMF THERAPY	
Intensity	Max. 24,700 Gs
Frequency	Max. 15 Hz

TECAR THERAPY	
Frequency	300 kHz & 448 kHz
Mode	Capacitive & Resistive
Output power	350 W

Detailed Specifications

Specification	Value
Power Supply	AC 100–240 V, 50/60 Hz
Output Power	10–3000 W (varies by module)
Frequency Range	1 Hz – 3000 Hz (varies by module)
Treatment Time	1–60 min adjustable
Display	10.1" Touch Screen
Dimensions (Single Unit)	450 × 360 × 160 mm
Weight (Per Box)	6–12 kg (varies by module)
Operating Temperature	10 °C – 40 °C
Storage Temperature	-10 °C – 60 °C
Certifications	CE, ISO 13485