

Low level steady heat compared to high level pulsed heat for the treatment of low back pain

Charles Chabal MD, PJ Dunbar MD, SA Stoner Ph.D. University of Washington, Evergreen Medical Center, Talaria Inc. Seattle, WA

Funding: NIH 2R44CA099305-02

Introduction

One third of all Americans suffer from back pain at some point during the year. Men and women are equally affected. The risk of experiencing low back pain from disc disease or spinal degeneration increases with age. Low back pain (LBP) is the most common cause of job-related disability and a leading contributor to missed work. Back pain is the second most common neurological ailment in the United States — only headache is more common (1). Heat has long been a mainstay treatment of low back pain. A number of recent well-designed studies demonstrated that heat reduces pain, improves function and may result in the use of fewer pain medications in sufferers of low back pain (2-5).

In spite of both empiric evidence and formal studies little is known about basic mechanisms surrounding heat and LBP relief. For instance, what temperature do subjects prefer, does the temperature response display dose curve characteristics, and how large of an area needs heated to produce an analgesic response? Based on earlier studies it was hypothesized that there was a dose dependent response curve associated with thermal analgesia (6). Furthermore the characteristics of the heating wave could additionally affect analgesic outcomes.

In a previous study we showed that subjects strongly preferred temperatures that were significantly hotter (mean 44.6 °C, range 42-48 °C) than those provided by chemical hot packs (39-40° C) and that heating a relatively small area of skin at high temperatures could produce analgesia comparable to that produced by heating large areas of skin at lower temperatures. In the present study we created an experimental device that allowed subjects to preselect levels of pulsed heat that were significantly higher than temperatures used in common commercial heating products. The hypothesis of this study was that the high-level pulsed heat delivered via the experimental device would be more effective than low-level continuous heat from a commercial hot pack in relieving LBP associated pain.

Procedures

This study compared the pain relieving effects of an experimental electronic heating device termed, *Series-3* with ThermaCare® wraps. The name Series-3 was used to designate the third laboratory iteration of the experimental device. Series-3 consisted of two electronic heating pads set to a temperature selected by the individual subject (figure 1 and 2). The temperature of the experimental device can be set by the user in a range from 42 to 47 °C via a wireless control. The Series-3 experimental device is shown in figures 1 and 2. The control device, ThermaCare®, is a commercially available, over the counter chemical hot pack. ThermaCare® is attached to the skin using its own elastic wrap and heats at a steady 40 °C.

All subjects met with a research assistant (RA) prior to the start of the study. The RA explained and demonstrated both devices, their purpose and the methods of the study. Subjects were randomly assigned to one of two groups, Series-3 or ThermaCare group. All subjects completed a questionnaire about their pain.

Once randomization occurred subjects rated their pain level using both a numeric pain scale and the lowa pain thermometer. Those subjects initially assigned to the control group had the device placed over their area of greatest pain. All ThermaCare wraps were allowed at least 30 minutes to warm up before being placed on the subject and a temperature of at least 40° was confirmed. Subject rated their pain levels at baseline and after 10 minutes, 20 minutes, and 30 minutes. After the first treatment session there was a 30-minute washout period after which subjects crossed over to the other treatment arm.

Subjects randomized to the experimental treatment group were shown the study device. The RA facilitated a run in period in which the subject was able to gradually increase the temperature of the Series-3 device starting at 108 degrees Fahrenheit (41 °C) up to a maximum of 116.6 degrees Fahrenheit (47 °C). Once the study temperature was selected subjects wore the device and provided pain assessments at baseline and after 10 minutes, 20 minutes, and 30 minutes.

After completing the study subjects filled out an exit interview questionnaire and were paid \$100 for study participation. Results were analyzed using paired t-tests.

Results

Pain Reduction: Heater Device and ThermaCare® device in LBP subjects

Results of a paired-sample t-test indicated significantly greater decrease in lowa Pain Thermometer scores from baseline to 30 minutes when participants used the Series-3 device ($t(43)$, $p<.001$). Similar differences were found from baseline to 10 minutes ($t(43)$, $p=.005$), and from the 20 to 30 minute assessment ($t(43)$, $p=.009$). No significant differences were found on reduction of scores on the IPT from 10 to 20 minutes.

Reduction of pain rating on the Numeric Rating Scale from baseline to 30 minutes was also greater with use of the Series-3 device ($t(43)$, $p<.001$). Those using the Series-3 device also reported greater reduction of pain on the Numeric Rating Scale from baseline to 10 minutes ($t(43)$, $p=.003$), and from 20 to 30 minutes ($t(43)$, $p=.003$). No significant differences were found from 10 to 20 minutes. See tables 1-4.

Conclusions

Both treatments produced significant reduction in pain in subjects who suffered from chronic low back pain. When compared to ThermaCare, the experimental condition produced faster onset and significantly better pain relief. In exit interviews subjects almost unanimously noted that they all preferred the warmer temperatures from the experimental device than that offered by the low level heat of the ThermaCare product. Many subjects also offered that they very much liked the pulsing sensation provided by the Series-3 device. Based on the rapid onset and heat profile it is likely that the analgesic response noted from the high-level pulsed heat may be related to stimulation of a-delta thermal receptors. The findings of this study may help identify unique properties and characteristics of a well-established treatment (heat) that can improve the speed of onset and quality of pain relief. Coupled with advances in wireless connectivity, lithium battery technology and low cost embedded chip sets these finding may offer novel mobile and more effective drug free options for patients who suffer from chronic pain.

References

1. National Institute of Neurological Disorders and Stroke, Low Back Pain Fact Sheet. www.ninds.nih.gov/disorders/backpain.
2. Akin, M. et al. (2001). Continuous low-level topical heat in the treatment of dysmenorrhea. *Obstetrics & Gynecology*, 97(3), 343-349.
3. Akin, M. et al. (2004). Continuous, low-level, topical heat wrap therapy as compared to acetaminophen for primary dysmenorrhea. *The Journal of reproductive medicine*, 49(9), 739.
4. Nadler, S. et al. (2003). Overnight use of continuous low-level heatwrap therapy for relief of low back pain. *Archives of physical medicine and rehabilitation*, 84(3), 335-342.
5. Nadler, S. et al. (2002). Continuous low-level heat wrap therapy provides more efficacy than ibuprofen and acetaminophen for acute low back pain. *Spine*, 27(10), 1012-1017.
6. Chabal C., Stoner S A., Dunbar PJ. (2013) Characteristics of Thermal Analgesia in Human Subjects *Pain Medicine* 14 (4), 557.

Timepoint	Series-3		ThermaCare	
	mean	median	mean	median
0-10 min	1.33	1.0	0.32	0.0
10-20 min	0.42	0.0	0.36	0.0
20-30 min	.70	1.00	0.54	0.0
0-30 min	2.44	2.00	1.21	1.0

Table 1: Mean and median change in participant pain LBP rating on a lowa Pain Scale at 10, 20, and 30 (N=27).

Timepoint	Series-3		ThermaCare	
	mean	median	mean	median
0-10 min	0.81	1.0	0.29	0.0
10-20 min	0.41	0.0	0.18	0.0
20-30 min	0.63	1.00	0.54	0.0
0-30 min	1.86	2.00	1.00	1.0

Table 2: Mean and median change in participant pain LBP rating on a Numeric Rating Scale at 10, 20, and 30 (N=27).

Timepoint	LBP N=27				
	Mean diff.	Std. Dev.	t	df	p
0-10 min	1.11	2.10	2.75	26	.01
10-20 min	0.4	1.26	.15	26	.88
20-30 min	.19	1.42	.68	26	.50
0-30 min	1.33	2.50	2.78	26	.01

Table 3. Mean differences: change in participant pain rating on lowa Pain Thermometer at 10, 20, and 30 minutes

Timepoint	LBP N=27				
	Mean diff.	Std. Dev.	t	df	p
0-10 min	.55	1.31	2.20	26	.04
10-20 min	.22	.58	2.00	26	.06
20-30 min	.11	.93	.62	26	.54
0-30 min	.89	1.43	3.24	26	.003

Table 4: Mean differences: change in participant pain rating on Numeric Rating Scale at 10, 20, and 30 minutes for participants with LBP using Series 3 and ThermaCare®.

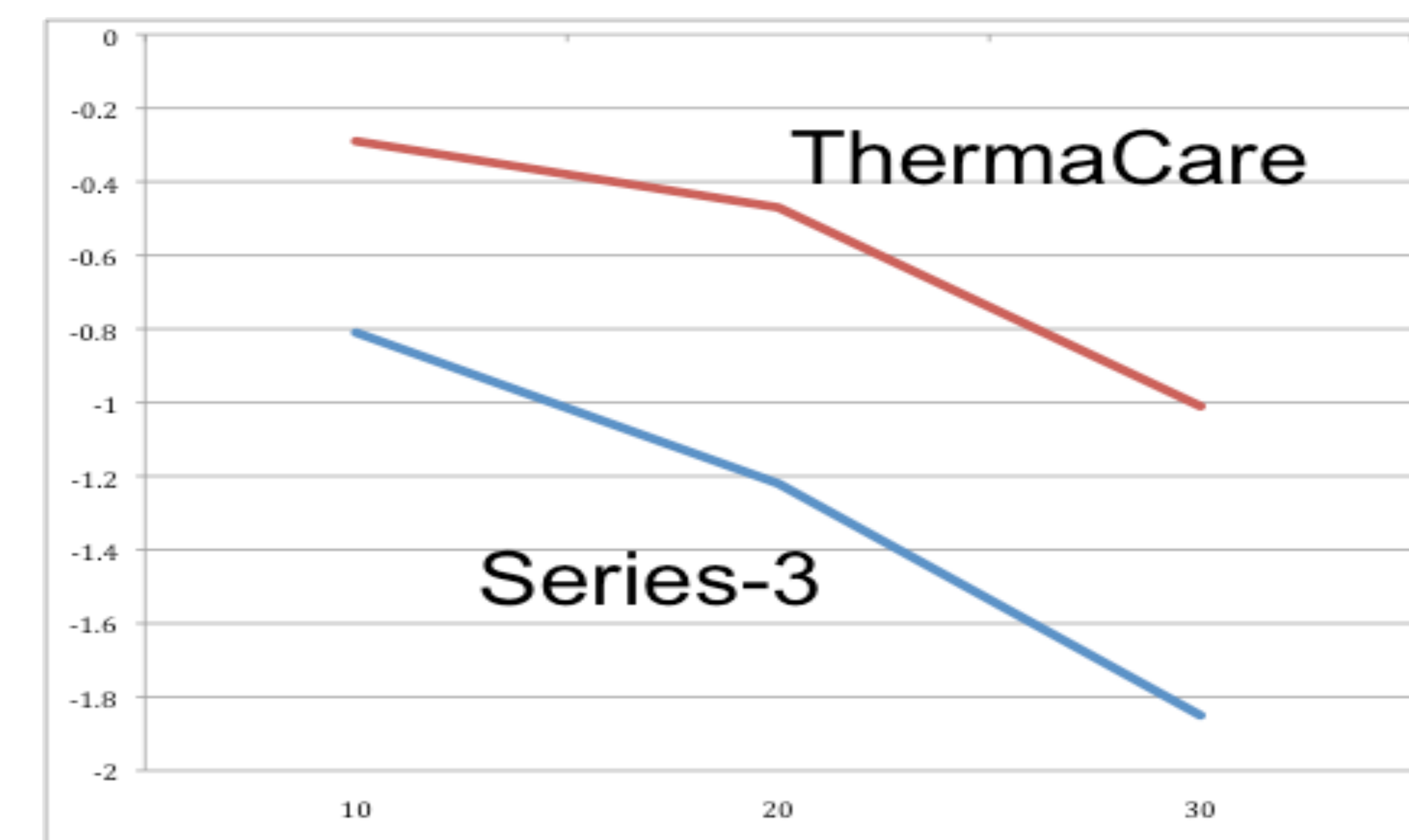


Figure 3: Reduction in LBP in subjects 0-10 scale after 10, 20, and 30 minutes of treatment.